
MT 410

IB MT 410 S-V10 2925, Revision Date: 2025-07-03

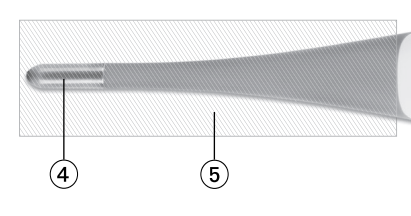Description of this thermometer

- ON/OFF button
- Display
- Battery compartment cover
- Measuring sensor / measuring tip
- Cleaning and disinfecting area (thermometer probe only)

This Digital Antimicrobial Medical Thermometer provides highly accurate readings over the human body temperature range. At the same time, it reduces the microbial flora and minimizes the dispersion of contagious microorganisms, providing high safety to the user.

Antimicrobial Copper Properties

Surfaces made or covered by special copper alloys, have strong antimicrobial properties against a wide variety of microorganisms*. Copper alloyed into the thermometer probe (see number 3) that while in contact with microbes and bacteria rupture their cellular membranes, thus destroying these microorganisms. This activity reduces the microbial flora on the coated area and - due to the «halo phenomenon» - simultaneously causes a drastic reduction in pathogens on the remaining body of the thermometer. Thermometers containing parts with antimicrobial copper alloys drastically reduce microbial flora, minimizing the dispersion of contagious microorganisms and thus providing high safety to the end user.

* *«Efstathiou A. Panos «The Role of Antimicrobial Copper Surfaces in Reducing Healthcare associated Infections», European Infectious Disease, Volume 5, Issue 2, Autumn 2011.*

Important safety instructions

- Follow instructions for use. This document provides important product operation and safety information regarding this device. Please read this document thoroughly before using the device and keep for future reference.
- This device is only to be used for measuring human body temperature through oral, rectal or axillary. Do not attempt to take temperatures at other sites, such as in the ear, as they may result in false readings and may lead to injury.
- Do not use this device if you think it is damaged or notice anything unusual.
- We recommend cleaning this device according to the cleaning instructions before first use for personal hygiene.
- The minimum measurement time until the beep is heard must be maintained without exception!
- Consider that different measurement locations may require continued measuring even after the beep, see section «Measuring methods / Normal body temperature».
- Do not attempt rectal measurements on persons with neural disorders. Doing so may aggravate or worsen the disorder.
- Ensure that children do not use this device unsupervised; some parts are small enough to be swallowed. Be aware of the risk of strangulation in case this device is supplied with cables or tubes.
- Do not use this device close to strong electromagnetic fields such as mobile telephones or radio installations. Keep a minimum distance of 3.3 m from such devices when using this device.
- Protect the device from impact and dropping
- Avoid bending the thermometer probe more than 45°!
- Avoid ambient temperatures above 60 °C. NEVER boil this device!
- Use only the commercial disinfectants listed in the section «Cleaning and disinfecting» to clean the device to avoid damage to the device.
- We recommend this device to be tested for accuracy every two years or after mechanical impact (e.g. being dropped). Please contact your MicroLife-Service to arrange the test.

WARNING: The measurement result given by this device is not a diagnosis! Do not rely on the measurement result only.

Batteries and electronic devices must be disposed of in accordance with the locally applicable regulations, not with domestic waste.

Read the instructions carefully before using this device.

Type BF applied part

Turning on the thermometer

To turn on the thermometer, press the ON/OFF button (1), a short beep signals «thermometer ON». A display test is performed. All segments should be displayed.

The last measurement reading will be shown on the display (2) automatically for 2 seconds. If the temperature is below 32 °C, the last temperature will be shown on the display (2). Then at an ambient temperature of less than 32 °C, an «LL» and a flashing «C» appear at the display field (2). The thermometer is now ready for use.

Function test

Correct functioning of the thermometer is tested automatically each time it is turned on. If a malfunction is detected (measurement inaccurate), this is indicated by «ERR» on the display, and a measurement becomes impossible. In this case, the thermometer must be replaced.

Using the thermometer

Before use, keep this thermometer out of any physical contact, for at least 2 hours. This utilizes the antimicrobial copper properties, as described in chapter «Antimicrobial Copper Properties». The use of antimicrobial copper is a supplement to and not a substitute for standard infectious control practices. Users must continue to follow all current infection control and cleaning practices. We recommend cleaning the thermometer as described in chapter «Cleaning and disinfecting».

Choose the preferred measuring method. When taking a measurement, the current temperature is continuously displayed and the «C» symbol flashes. If the beep is heard 10 times and the «C» is no longer flashing, this means that the measured increase in temperature is less than 0.1 °C in 5 seconds.

To prolong the battery life, turn off the thermometer by briefly pressing the ON/OFF button (1). Otherwise the thermometer will automatically turn off after about 10 minutes.

Storage of Measured Values

If the ON/OFF button (1) is pressed for more than 3 seconds when turning on the thermometer, the automatically stored maximum temperature during the last measurement will be displayed. At the same time, a «M» for memory will appear on the display. About 3 seconds after the button is released, the temperature value disappears and the thermometer is ready for measurement.

Measuring methods / Normal body temperature

► In the armpit (axillary) / 34,7 - 37,3 °C

Wipe the underarm with a dry towel. Place the measuring sensor (4) under the arm into the center of the armpit so the Measuring sensor is touching the skin and position the patient's arm next to the patient's body. This ensures that the room air does not affect the reading. Because the axillary takes more time to reach its stable temperature wait at least 5 minutes, regardless of the beep sound.

► In the mouth (oral) / 35,3 - 37,5 °C

Do not eat or drink anything hot or cold 10 minutes before the measurement. The mouth should remain closed up to 2 minutes before starting a reading.

Position the thermometer in one of the two pockets under the tongue, to the left or right of the root of the tongue. The measuring sensor (4) must be in good contact with the tissue. Close your mouth and breathe evenly through the nose to prevent the measurement from being influenced by inhaled/exhaled air.

If it is not possible due to blocked nostrils, another method for measuring should be used.

Aprox. measuring time: 30 seconds

► **In the anus (rectal) / 36,6 - 38,0 °C**
The use of a probe cover and the use of a lubricant is recommended. Carefully insert the measuring sensor (4) of the thermometer 1 to 2 cm into the anal aperture.

If you are unsure of this measurement procedure, you should consult a professional for further information.

Aprox. measuring time: 30 seconds

Cleaning and disinfecting

For disinfection in home use environment, use a 70% Isopropyl alcohol swab, or a cotton tissue moistened with 70% Isopropyl alcohol to wipe surface pollutants off the thermometer probe (note: consider the application and safety instruction of the disinfectant manufacturer). Always start wiping from the end of the thermometer probe (approx. at

microlife

the middle of the thermometer) towards the thermometer tip. Afterwards the entire thermometer probe (see number 3) in the (drying) should be immersed in 70% Isopropyl alcohol for at least 5 minutes (max. 24 hours). After immersion, let the disinfectant dry off for 1 minute before next use. Avoid immersing or wiping the display to protect it from fading.

The thermometer is not intended for professional use.

Battery installation and replacement

When the «▼» symbol (upside-down triangle) appears in the display, the battery is fast and needs replacing. To replace the battery remove the battery compartment cover (3) from the thermometer. Insert the new battery with the «+» at the top. Make sure you have a battery of the same type and brand. Batteries can be purchased at any electrical store.

Technical specifications

Type:	Maximum thermometer
Measurement range:	32.0 °C to 42.9 °C
Temp. < 32.0 °C: display «LL» for low (too low)	
Temp. > 42.9 °C: display «H» for high (too high)	
± 0.1 °C; 32.0 - 33.9 °C; 42 °C	
± 0.2 °C; 34.0 - 33.9 °C and 42.1 - 42.9 °C	
10 - 40 °C; 15 - 95 % relative maximum humidity	
Storage conditions:	-25 - +60 °C; 15 - 95 % relative maximum humidity
Battery:	SR41 (1.55V)
Expected lifetime:	approx. 2700 measurements (using a new battery)
IP Class:	IP22
Reference to standards:	EN 12470-3, clinical thermometers; ASTM E1112; IEC 60601-1; IEC 60601-1-2 (EMC); IEC 60601-1-11
Expected service life:	5 years or 10000 measures

This device complies with the requirements of the Medical Device Directive 93/42/EEC. Technical alterations reserved. **Guarantee** This device is covered by a 5 year guarantee from the date of purchase. During this guarantee period, at our discretion, MicroLife will repair or replace the defective product free of charge. Opening or altering the device invalidates the guarantee. The following items are excluded from the guarantee:
• Transport costs and risks of transport.
• Damage caused by incorrect application or non-compliance with the instructions for use.
• Damage caused by leaking batteries.
• Damage caused by accident or misuse.
• Packaging/storage material and instructions for use.
• Regular checks and maintenance (calibration).
• Accessories and wearing parts: Battery.
Should guarantee service be required, please contact the dealer from where the product was purchased, or your local MicroLife service. You may contact your local MicroLife service through our website: www.microlife.com/support. Compensation is limited to the value of the product. The guarantee will be granted if the complete product is returned with the original invoice. Repair or replacement within guarantee does not prolong or renew the guarantee period. The legal claims and rights of consumers are not limited by this guarantee.

This device is covered by a 5 year guarantee from the date of purchase. During this guarantee period, at our discretion, MicroLife will repair or replace the defective product free of charge. Opening or altering the device invalidates the guarantee. The following items are excluded from the guarantee:
• Transport costs and risks of transport.
• Damage caused by incorrect application or non-compliance with the instructions for use.
• Damage caused by leaking batteries.
• Damage caused by accident or misuse.
• Packaging/storage material and instructions for use.
• Regular checks and maintenance (calibration).
• Accessories and wearing parts: Battery.
Should guarantee service be required, please contact the dealer from where the product was purchased, or your local MicroLife service. You may contact your local MicroLife service through our website: www.microlife.com/support. Compensation is limited to the value of the product. The guarantee will be granted if the complete product is returned with the original invoice. Repair or replacement within guarantee does not prolong or renew the guarantee period. The legal claims and rights of consumers are not limited by this guarantee.

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